




MEDCONCEPT

BDS-05

BDS-05 Mobile Cart Veterinary Dental Unit is an excellent workstation for pet animal dental cutting, grinding, scaling, polishing, root canal treatment, etc. Providing 5 integrated handles that cover the main functions required by veterinarians, BDS-05 can completely fulfill the needs of animal dental treatment, allowing veterinarians to offer the most comprehensive and highest quality dental care.

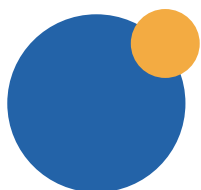

Product Features

Smart

Small area size 40*43cm, saving space:
Adjustable height 90cm - 105cm, suitable for the small space clinics.

High Efficiency & Low Noise

High-quality low noise oil-free air compressor: Water flow pressure can be simply and efficiently controlled

LED Lighted Ultrasonic Scaler

Plug-in handle with LED light for clearer operation vision

3-way air/water Syringe

360° swivel handpiece, needle can be autoclaved

Low-Speed Handpiece


Low-speed micro motor with external water channel

High-Speed Handpiece

Stable speed with water column cooling

Suction Handpiece

Easily suck and clean the liquid in mouth

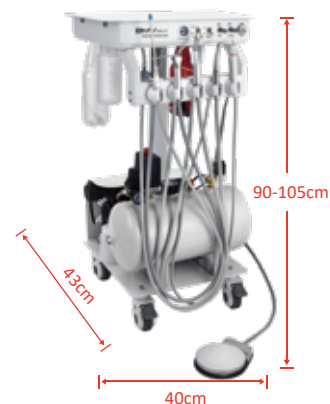

Technical Parameters

Items	Specification
Compressor	High-quality low noise oil-free air compressor
Input Power	2.6-2.9A, 220/230V, 50/60Hz
Output Pressure	0.45~0.5MPa
Gas Tank	12L
Sound Level	<65dB
High-Speed Handpiece	330000-390000 rpm
Low-Speed Handpiece	5000-20000 rpm
Ultrasonic Scaler	LED lighted, Tip ultrasonic frequency: 28kHz+3kHz
Dimension	43 *40 *95 Height adjustable from 90cm -105cm
Weight	28kg

MEDCONCEPT

MEDCONCEPT TECHNOLOGY LIMITED

Address Shenzhen:

B301, Languang Technology Park, No. 27, Gaoxin North Sixth Road, Nanshan District, 518057, Shenzhen

Tel/Fax: +86-755-23732826

Mail: maxim@medconcept-tech.com

Web: www.medconcept-tech.com